

LED Beauty Mask

User Manual



Catalogue

1.	Product principle/productfunction	P - 01
2.	Product parts introduction	P - 03
3.	Use method/cautions	P - 04
4.	Maintenance and storage/people unsuitable to use	P - 05
5.	Product specifications/troubles shooting	P - 06
6.	Warranty card	P - 08

Please read this manual before use

Thank you for purchasing the spectrum rejuvenation instrument. Please read the operation manual carefully and pay special attention to safety precautions before use, in order to operate and use the product correctly.

This manual provides users with safe and effective ways of operation. Please keep it properly for future reference after reading.

Product principle


This product selects the red, blue, orange light source, through the nanotechnology light guide material, uses the biological phototherapy principle, to the skin gentle massage, lets the skin metabolism, causes the skin to be transparent and luster, it can effectively reduce acne, stimulate collagen production, improve cell microcirculation and combine the essence of facial mask to achieve the effect of skin beautification.

Product features

LED photon therapy is widely used by skin specialists because it's a noninvasive painless skin treatment with significant effect and no chemical additives or pollution. The wavelength of 620-630nm released by the red light can promote blood circulation, accelerate metabolism, and promote the growth of subcutaneous collagen, so as to eliminate wrinkles and revitalize the skin. The wavelength of 460-470nm released by the blue light can effectively eliminate acne bacteria so that the acne can be eliminated from the root, improve the environment inside the pores and restrain the bacteria deep in the skin. The 580-590nm wavelength released by the orange light can improve the alternating function of cell oxygen, improve microcirculation, dilute pigment, reduce skin roughness and poor circulation.

1. Promote skin to produce collagen, make skin more bright and translucent, reduce fine lines and wrinkles.
2. Tighten pores, improve skin brightness and make skin glossy
3. Dilute freckles and red blood cells, accelerate the improvement of skin and achieve better skin repair after sunburn.
4. Provide skin cells with activation energy, accelerate the metabolism of the skin, and keep the skin as supple as baby.

Product Parts Introduction

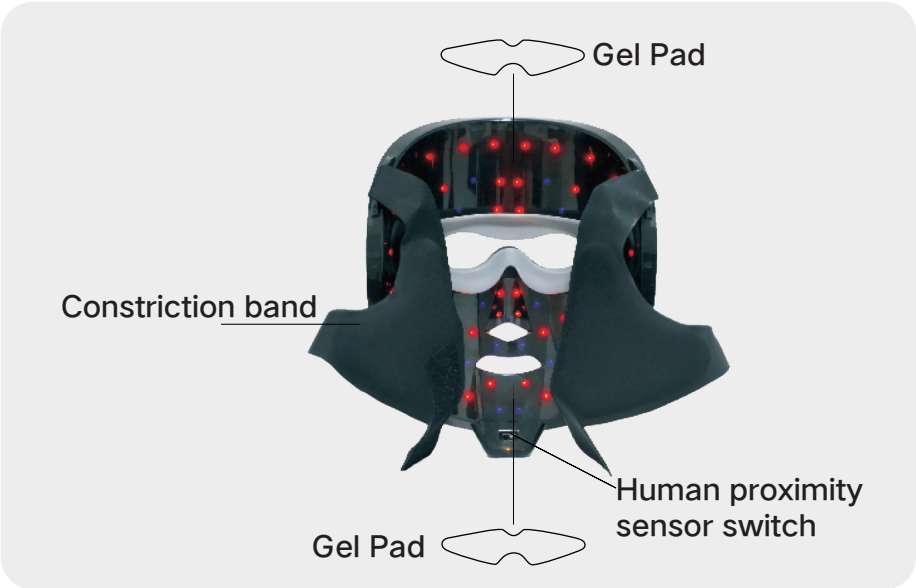

Machine Condition	Charging state	Battery status	Indicator status
Start up	Uncharged state	Low Battery	The white light are flashing
Start up	Uncharged state	High Charge	The white light is always on
Start up	Charging state	Be Full of	The white light is always on
Start up	Charging state	Charging process	White breath light
Shut down	Charging state	Charging process	White breath light
Shut down	Charging state	Be Full of	White light out
Shut down	Uncharged	/	White light out

- 1.Please stay away from fire or corrosive liquids.
- 2.Please do not disassemble and repair this product by yourself.
- 3.Keep out of reach of children.
- 4.Do not apply heavy pressure to avoid deformation and damage.

Maintenance and storage

- 1.After use, wipe the host machine with a dry cloth,do not wash it with water directly.
- 2.If this product is damaged, please do not continue to use.
- 3.Disassembly of this product is strictly prohibited.In case .
- 4.Store in a cool, dry place, away from moisture and strong sunlight .
- 5.Please do not put the product into the liquid.
- 6.Keep away from children in case of fear of danger.

Use of taboo groups

- 1.It is forbidden to use in pregnant women, and the Nano light source will not affect the fetus, but it is still not recommended.
- 2.People with epilepsy are forbidden to use it .
- 3.It is forbidden to use in patients with thyroid diseases and thyroid drug users.
- 4.Not For people with light allergies.

Product Specification

Product name	S520
Number of LED lights	66 Only
Product size	260X208X100mm
Battery capacity	2000mAh
Net weight of product	398g
Product material	Silica Gel+ABS+PVC
Electric Current	Less than or equal to 1.8A
rated voltage/Charging time	5V/About three hours
Rated power	Less than or equal to 6W
Optical wavelength	Red: 620-630NM; Blue: 460-470NM; Orange: 580-590NM
Applicable skin	All Skin
Accessories	Power Cable; manual

Trouble shooting

Malfunctions	Cause of failure	Fault handling
No work	It's not on	Long press to touch the switch to open
	Low Battery	Check the charger and look atthe Indicator light for charging
	Bad socket contact	Please reinsert or replace the charger
	Automatic shutdown	Press the touch switch toturn it on again

Benefits of Uno Vita mask

The far-infrared wavelength of 6~14 μ m generated when graphene mask heated is quite similar to the human body and can achieve the same frequency resonance effect with human body, thereby fundamentally restoring the skin barrier health, better absorbing skincare nutrition and improve skin problems.

The LED high-tech photon skin rejuvenation beauty instrument uses the principle of photodynamics to activate deep cells and make the skin better metabolize. After the light is absorbed by the skin, the light energy is converted into intracellular energy, which relaxes and strengthens the capillaries and produces a photochemical reaction-enzymatic reaction with the skin, increasing the content of glycogen and protein, thereby promoting cell metabolism and synthesis

Warrantycard

Product Name		LED Light Therapy Mask	Model	
Distributor information				
Buyer information	Name:			
	Address:			
	Tel:			
Cause offailure				

* Please fill out the warrantycard in detail when repairing. Thankyou for your!

