# **User Manual**



## **Red Therapy Mask**

## **User Guider**

Create a better life with LED technology

## Content

| Product Principle                     | . 2 |
|---------------------------------------|-----|
| Instructions                          | . 3 |
| Product Details                       | . 4 |
| Parameters                            | 5   |
| Mode Optional7                        | '-8 |
| How to recharge for the mask and neck | 9   |
| Precautions & Forbidden to Use Crowd  | 10  |

### NOTE

Thank you for purchasing the Silicone LED Mask. Please read the guidebook carefully before using it, and pay attention to the Precautions in order to operate and use this product correctly.

This manual provides you with safe and effective operation methods. Please keep it properly after reading it for future reference.

#### **Product Principle**

WMA has confirmed that the use of specific wavelengths of light can enhance the vitality of skin cells and tissues through close irradiation, and better improve the skin condition. the Silicone LED Mask uses the principle of biophototherapy with nanotechnology to guide light materials, won't generate high heat, and uses low-energy pure light and touchable light adjustment technology to care for the skin. The treatment has the characteristics of natural mildness, remarkable curative effect, low risk, no side effects, etc. It is a skin care expert and is suitable for all types of skin.

The LED silicone color light skin rejuvenation mask uses the principle of photodynamics to activate deep cells and allow the skin to better metabolize. After the light is absorbed by the skin, the light energy is converted into intracellular energy, which relaxes and strengthens the capillaries and produces a photochemical reaction-enzymatic reaction with the skin, increasing the catalase that acts as the "guard" and "cleaner" of cells and the activity of SOD. Increase the decomposition of ATP, which is the source of cellular energy.

Promote cell metabolism and synthesis, increase glycogen and protein content, and stimulate 5 times collagen. Stimulate the secretion of BFGF and EGF in the body, rearrange collagen fibers and elastic fibers, block and eliminate the formation of melanin, accelerate cell growth, accelerate blood circulation, and stimulate fibrous cell tissue It produces collagen, increases skin elasticity, and plays an important role in repairing aging skin, skin with sores, diminishing spots, lifting firmness, and relieving sunburns.

#### Instructions

- 1. Clean the skin of the face and dry the water stains on the face with a soft towel
- 2. Take out the product and connect the usb line with the controller
- 3. Turn on the Power Button

4. Press Color Button choose the mode, then press the Time Button to adjust the working time as you want, 15–30 min is the best if need recommend. And don't stare at the light for too long

5. Short Press the Power Button can adjust the strength of the light


6. It is recommended to use this product while lying on the bed or sofa to relax the whole body

7. It will have a better results if use with facial mask and essence or lotion

8. After use, take off the silicone skin rejuvenation mask and do a good job of cleaning, so that you can continue to use it next time.



#### **Product Details**

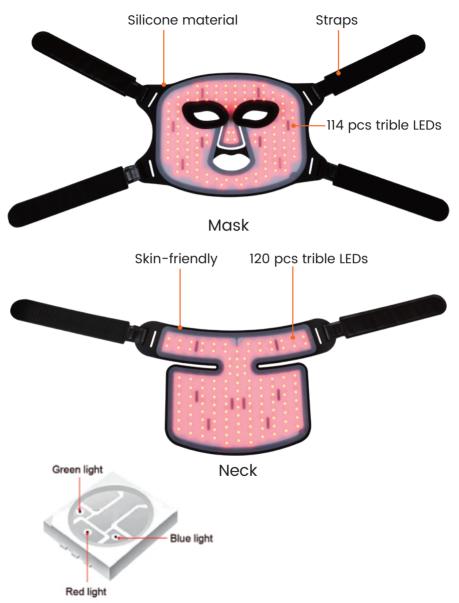

#### Parameters

| Product Name:        | Skin LED Light Therapy Face Mask      |
|----------------------|---------------------------------------|
| LED Qty:             | 234pcs                                |
| Product Size:        | 320*210mm (Mask) 、 330*210mm (Neck)   |
| Wavelength:          | 470nm 630nm 520nm                     |
| Brightness Function: | 0%-100%                               |
| Color:               | Black                                 |
| Net Weight:          | 680g                                  |
| Mode:                | 7+1 modes                             |
| Power:               | Max 12W                               |
| Material:            | Sillicone                             |
| Timer:               | 0-30mins, quickly double-click to set |
| Frequency:           | 0-10000Hz                             |
| Application:         | For face and neck                     |
| Customize:           | Logo/shape/wavelength/color/box       |
| Warranty:            | 1 year                                |



1\*Light For Mask





1\*User Manual

3\*Fur Straps



1\*Light For Neck







3\*Hook Straps

#### **Operating Instructions**

 Display Window: Shows all the function of interface(Timer, Mode, Brightness, Pulse function).

2. SET Button: Setting timer, wavelength selection , brightness and Pulse



3. 4: "+""-" Button : Press to set the time duration
from 1-30minutes, in 1 minute increments, the default
timer is 10minutes.
Press to set the brightness of light, allowed you set it
from 1% to 100% of light intensity.

Press to choose different 8 modes: red, green, blue, red & green, red & blue, green & blue, red & green & blue, loop(7 types of light cycle on)

Press to set the pulse function. 1 - 9999Hz, NULL

#### 5.Turn on

When setting well the Timer, Mode, Brightness, Pulse at above steps, then turn on the button.

#### **Mode Optional**

**Red:** It can improve cell activity, speed up metabolism and promote collagen formation. It has the effects of whitening and lightening, rejuvenating and wrinkling, increasing skin elasticity, improving skin condition, repairing damaged skin, reducing pores and tightening skin, anti-aging and anti-oxidation.







Green: It can reduce the secretion of skin oil, balance the skin, effectively relieve stress, dredge lymph and drain swelling.







Blue: It has the effect of bactericidal, anti-inflammatory and inhibiting inflammation, and can remove acne and sores, sterilize and reduce inflammation, repair skin, and leave no scars and pigments.







Red & Green: Replenish energy to skin cells, decompose pigments, promote lymphatic drainage, improve skin roughness and wrinkles, treat skin redness, fever, ringworm, etc., and enhance immunity.







Red & Blue: It is a dual-frequency light of red light and blue light. It combines the effects of two kinds of phototherapy, especially in the treatment of acne and acne scars and has a particularly good effect and repairing effect.







Green & Blue: It can increase cell energy step by step, and has a good effect on metabolism.



Red & Green & Blue: Penetrates the deep layers of the skin, accelerates skin metabolism, breaks down pigment spots, improves fine lines and sagging skin.





Loop (7 types of light cycle on): All benefits




#### How to recharge for the mask and neck



2

#### Power supply for power bank

Connect the charging cable from the mask and neck part to the power bank.

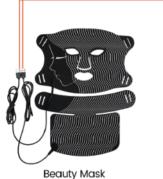



Power Bank&Controller

#### Charging the controller

20 ſĽ. × / Power Bank&Controller Charging Cable Plug Socket

Use Type-C to charge the controller.

Notice: After charging for 3 hours, it can be fully charged and can be used for 3.5 hours.

#### **Precautions**

1. Try to avoid looking directly at the light source part of the mask or staring at the light source for a long time. Keep your eyes closed during use. If you feel the light is too strong, you can cover it with cotton pads or eye masks.

2. It is recommended to use it with natural disposable facial masks with good quality ingredients or skin care products such as acne and whitening, which will help the absorption effect.

3. Because the mask is made of silicone plastic, there may be an odor after opening the package, which is a normal phenomenon, and the odor will disappear soon.

4. There are circuits and LED luminous bodies inside the luminescent film, which cannot be squeezed, folded, or pressed.

5. When using in bed before going to bed, it is best to put it back in the box after the treatment is used up, or at least take it off the bed.

6. When cleaning the luminous mask, it is strictly forbidden to soak in water, just wipe it with a damp cloth, and the elastic band can be washed with water.

7. Do not scrub the mask with corrosive chemical solvents.

8. Avoid exposure to sunlight.

#### Forbidden to Use Crowd

- $\cdot$  Do not use for pregnant women, children, sheep with epilepsy;
- · Do not use for thyroid patients or thyroid medications;
- Do not use for photo allergies or taking drugs that can cause photoallergies;
- · Do not do for those with skin trauma;

• Normal People who are uncomfortable should stop using it. After asking the doctor, follow the doctor's advice.



Our products are not intended to diagnose, treat, cure, or prevent any disease. Our products are low-risk, general wellness / fitness products in accordance with the "General Wellness: Policy on Low Risk Devices" draft released January 20, 2015.